





- excellent outdoor performance
- manoeuvrable and compact

SIGMA - midwheeldrive

The SIGMA gives you the freedom to go anywhere you want. It combines all the best features in one electronic wheelchair.

Thanks to the SIGMA's unique suspension and the LiNX steering, driving becomes an amazingly smooth experience and delivers an excellent performance in outdoor usage.

Even when fully equipped, the SIGMA remains compact and highly manoeuvrable. Perfect for negotiating sharp turns in tight situations as well as in open spaces.



The SIGMA wheelchair is available in different configurations and is built to measure. We strive to give you the best possible experience and the ultimate seating and driving comfort. A wide range of options allows the SIGMA to be adjusted to your personal needs and lifestyle.

- unique suspension
- driving comfort

SIGMA - VERMEIREN design

Our professional designers developed a unique suspension system with 3 pendulum arms. As such, all 6 wheels are in constant contact with the ground... on all sorts of terrains.

Because of this unique 3-points suspension, the SIGMA ensures stability at any given time and can handle height differences up to 50 mm (65 mm for the SIGMA 230).

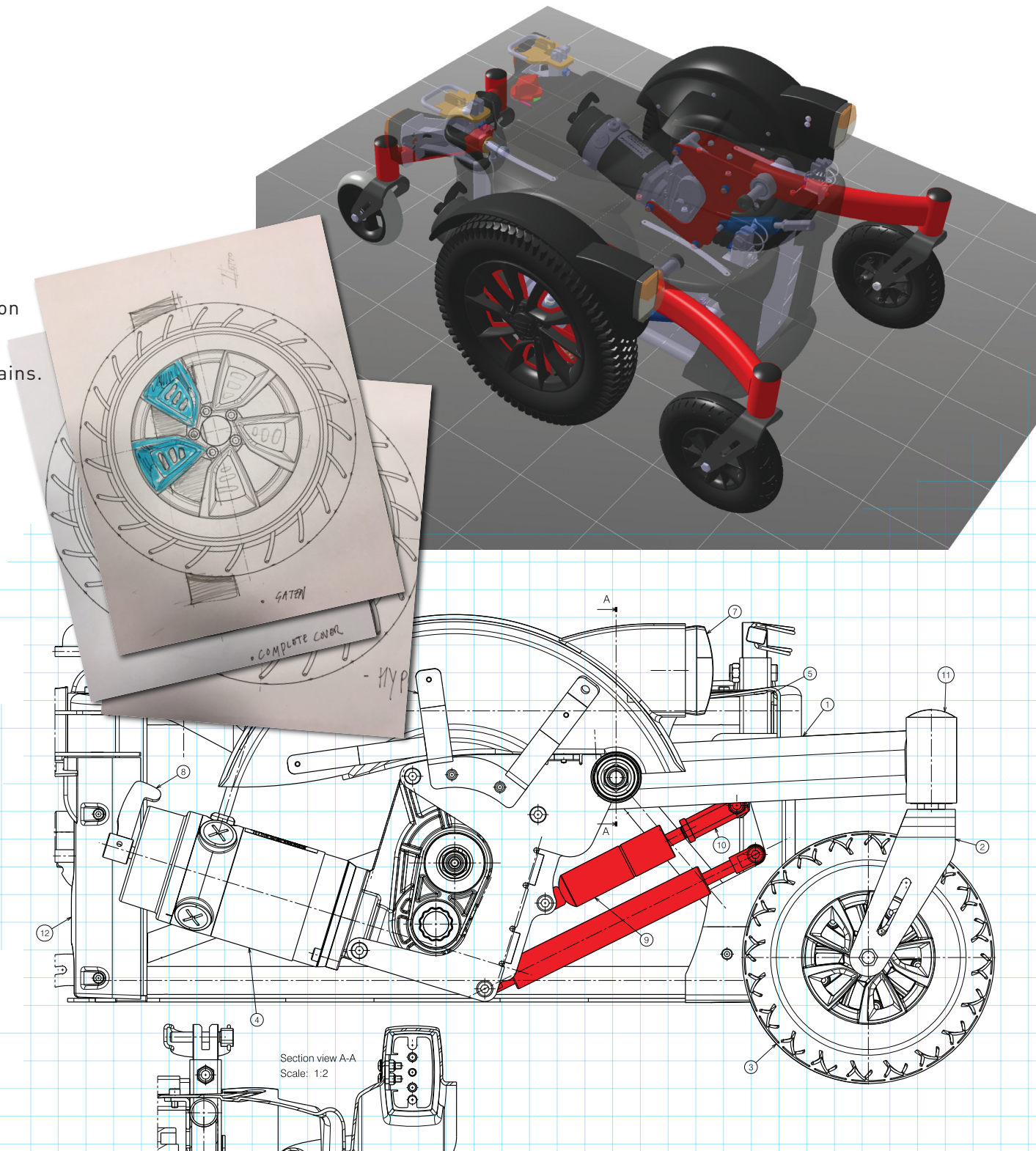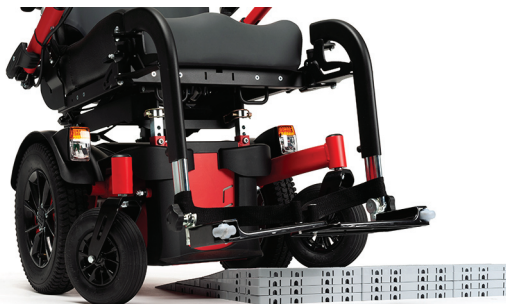



The next generation wheelchair control system

- A seamless and expandable system
- Modular components – mix and match modules
- State-of-the-art technology
- Revolutionary approach to programming



 adaptive controllers

SIGMA - LiNX controllers

The SIGMA comes with a **LiNX** controller system.

LiNX is an adaptive steering which constantly adjusts itself towards the driving conditions. It will detect when the SIGMA is driving on two different undergrounds (for example grass and concrete) and will compensate by adding extra power to one side as to keep the balance of power even and keep the SIGMA going in a straight line.

Due to its modular construction, the **LiNX** controller can easily be upgraded by adding options, adjust some of its functions and conduct maintenance.

One of the most used controllers is the basic controller, a standard controller with access to most of the SIGMA's features.




SE08 - Extra controller




SE09 - Attendant controller




SE64 - 10-way switch

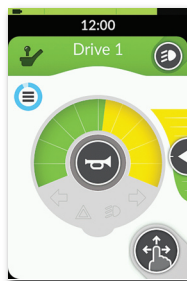


LiNX controller with LCD touch screen
Allows access to a large number of functions (inclination of seat and backrest, lighting ...) with direct visualization on the LCD screen for more comfort and practicality.

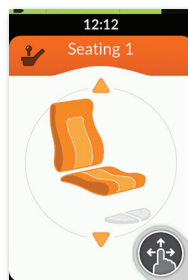


Seat adjustments

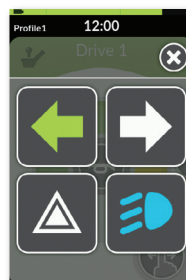


Visible and clear icons

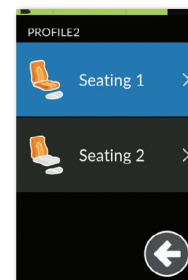


Personalised seating

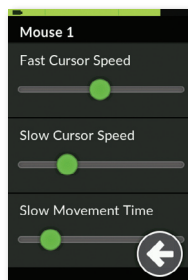

Fine-tuning of the settings

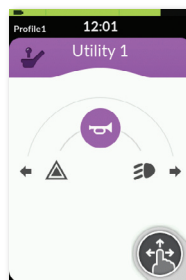

Different utility programs

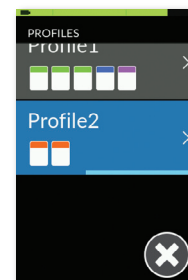

Personalised profiles

The **LiNX** interface is available on iOS and PC, allowing you to configure a chair or analyse its data quickly and easily.

LiNX offers comprehensive data and statistics to monitor the performance of your wheelchair. They can be viewed in real time with a Bluetooth® connection.

The dongle key can be used to adjust any chair equipped with a **LiNX** controller.



SE83 - Mini joystick



SE82 - Button controller

SIGMA - LiNX additional controllers

Additional controllers consist out of a separate screen, an extra controller that can be mounted at the back, 10-way switch for fingertip control and a number of special adaptive controllers to find the most ideal solution for each user.



■ freedom of movement

SIGMA - controllers

VERMEIREN strives to enhance the personal comfort of each users and offers different types of controllers and systems depending on the personal needs of each person. From basic controllers to made to measure chin controllers with fully integrated touchscreens.

SE04 - Chin controller



SE63 - Softball for joystick



SE02 - LiNX basic controller



SE02 - LiNX LCD controller



SE53 - Swing away controller



SE53A - Swing away controller
with lock positions



SE53Z - Swing away down controller



SE55 - Reinforced controller

personal comfort

SIGMA - Seating comfort

Comfortable backrests are essential to any persons daily comfort. From soft, preshaped backrests to backrests with special ventilation pores or the more stylish, artificial leather ones... all types are compatible with the SIGMA wheelchair and will add to a more pleasant experience when driving it.



Armrest in PU material



Armrest in artificial leather



Armrest in soft material



Back and seat in perspiration ventilating material



Backrest and seat in SOFT material



Backrest and seat in artificial leather

■ seating configurations

Optional seating configurations are available upon request. From an optional lift (SE42) to bridge a gap up to 300 mm or additional seat tilting up to 50° (SE43). These extra functions add to the quality of your lifestyle and can make things easier.



■ modern design



SIGMA 230

Design does not always mean you have to give in when it comes to functionality. For the users that demand a little extra, we have created the SIGMA 230. Equipped with bigger front wheels, steel front forks and an additional obstacle height up to 65mm.





SIGMA - Your guarantee of quality

Regarding quality, the SIGMA complies with the highest norms, laid down by certified test institutes in most European countries.

The wheelchair must be adapted to you, not vice versa

Our basic principle is: the wheelchair must be adapted to you, not the other way around. That is why our wheelchairs are evolutive and an answer to your specific needs.

All types of wheelchairs can be adapted, have numerous options and can be fully personalised so that there is a wheelchair for every person. This is the result of our continuous research and experience since 1957.



options

B06 - manual legrests

BZ8 ART - manual comfort legrests

BZ8E ART - electrical comfort legrests

B15A - table in ABS material

B31 - crutch holder

L04 - side support

SE24 - cup holder

SE26 - push bar

BZ10E - electrical footplate

L55 - adjustable headrest

L58 - adjustable comfort headrest

L44 - tension adjustable backrest

B58 - safety belt with AUTO lock

SE16 - plexiglass side plate

SE27 - topcase

SE74 - side mirror



Block profile



Block profile



Line profile



Puncture free





B06



BZ8 ART



BZ8E ART



B15A



B31



L04



SE24



SE26



BZ10E



L55



L58



L44



B58



SE16



SE27




SE74



SIGMA - matte and glossy colours

SIGMA is available in a wide range of different and vibrant colours to match your personal style. Either go for a subtle matte colour or a more striking glossy colour. Your colour of choice is represented in all the details of the SIGMA including the inner side of the central wheel rim and the lower front side of the battery encasing.

matte colours

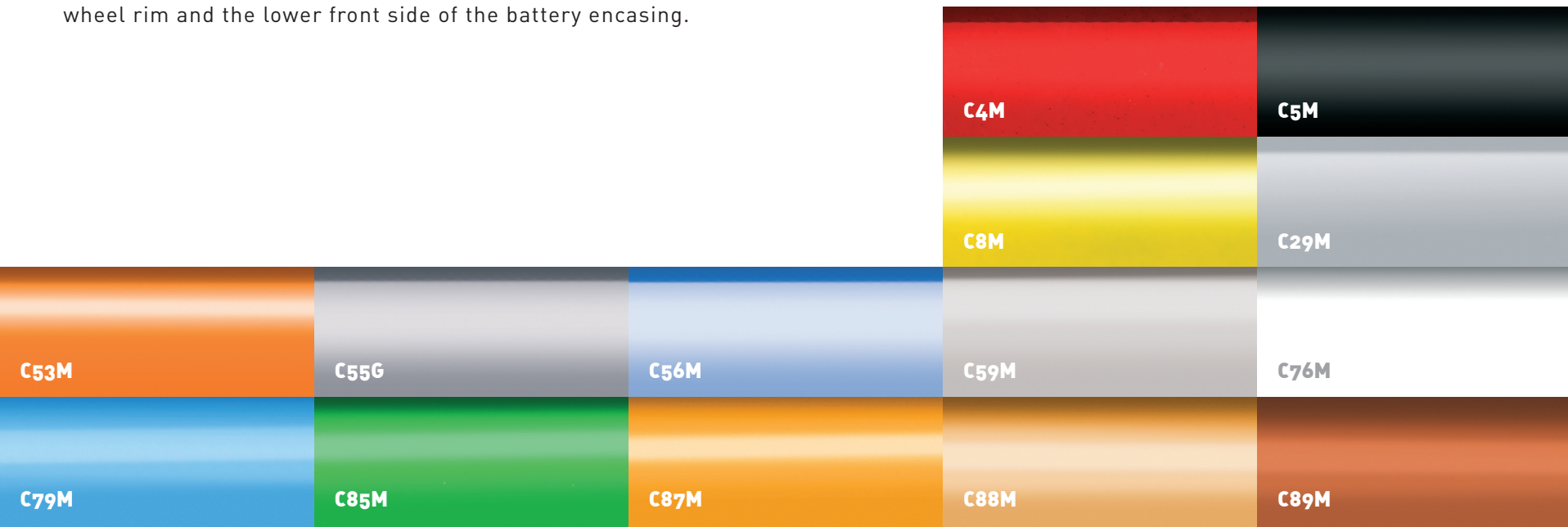

dimensions

Length, width, height... how much does it weigh? Do I need a bigger or a smaller seat width?... All these are important questions that we answer in our technical manuals to provide everyone with accurate information that will help you to make the right decisions on which model to get.

Description	Dimensions per configuration with BZ10E and L34 TUB			
	SE30	SE39	SE43	SE42
Max. weight of user	140 kg			
Wheelchair class	B			
Overall length with leg rest	1177 mm		1163 mm	1198 mm
Overall width	Sigma Sigma 230 618 mm - 618 mm - 630 mm - 660 mm - 680 mm 672 mm - 672 mm - 672 mm - 672 mm - 680 mm			
Overall height	1055 - 1115 mm		1075 - 1135 mm	1100 - 1400 mm
Total mass	141 kg		154 kg	168 kg
Mass of the heaviest part	133.75 kg		146.75 kg	160.75 kg
Masses of removable parts	Footrest B06 Backrest Charger	1,95 kg (per piece) 3.35 kg 0,95 kg		
Static stability downhill		6°		
Static stability uphill		6°		
Static stability sideways		6°		
Energy consumption *		Drive A: 35 km Drive B: 35 km		
Dynamic stability uphill		6°		
Obstacle climbing	Sigma Sigma 230	50 mm (with a run-up of 0.5m) 65 mm (with a run-up of 0.5m)		
Maximum speed forward		Drive A: 6 km/h Drive B: 10 km/h		
Max. braking distance from max. speed		6km/h: 1 m , 10 km/h: 2,1 m		
Seat plane angle	0 - 21° (manual)	0 - 20° (electric)	0° - 50° (electric)	0° - 45° (electric)
Effective seat depth	420-440-460-480-500-520 mm			
Effective seat width	400-420-450-480-500 mm			
Seat surface height at front edge (adjustable)	435-450-465-480-495 mm		453 to 513 mm	478 to 778 mm
Backrest angle	0° - 30°			

Description	Dimensions per configuration with BZ10E and L34 TUB			
	SE30	SE39	SE43	SE42
Backrest height	610 mm			
Footrest to seat distance (BZ10E)	365 - 465 mm			
Leg to seat surface angle	90° - 153°			
Armrest to seat distance	250 - 530 mm			
Minimum turning radius	1200 mm			
Ground clearance (bottom base frame)	65 mm			
Diameter mid wheels	350 x70 mm			
Type of mid wheels	3.00-8, air tyre			
Diameter front wheels	Sigma Sigma 230 200 x 50 mm 230 x 70 mm			
Type of front wheels	Rubber tyre			
Diameter balance wheels	140 x 40 mm			
Type of balance wheels	Rubber tyre			
Base frame with wheels – dimensions	Sigma Sigma 230 962 x 618 x 377 mm 962 x 672 x 377 mm			
Description electrical	Drive A		Drive B	
Battery (Non-spillable AGM battery, maintenance-free)	Capacity: Output:		80 Ah (C20) 12V DC	
Driving motors	2 x 220W with EM brake		2 x 350W with EM brake	
Battery charger	Input: Max. charging current: Protection class Insulation class II		220-240V AC, 50/60Hz 6 A IP21 II	
Power module	Blow fuse		150A	

* Theoretical range, calculated in accordance with EU standards under ideal conditions. Speeds and range may deviate (approx. 15%).

Measurements in mm and kg. Deviations in sizes and weight (15 mm and 1,5 kg) are possible. Subject to technical alterations. Images for illustration purposes only. Certain features may be optional.



- unique suspension
- driving comfort
- excellent outdoor performance
- manoeuvrable and compact
- extensive product range

distributed by

