

Forest 3

ADULT POWERED
WHEELCHAIR



VERMEIREN
we care for you

Forest 3

DRIVING SAFELY AND COMFORTABLY

Known for its axle system, the Forest 3 allows its user to drive safely and comfortably on various surfaces. With Forest 3, you'll overcome obstacles with ease in your daily life. A wheelchair made for adventures in- and outdoors.



COMPACT FOR EASY MANEUVERING

The small and compact size of Forest 3 gives you a smooth turning radius. Perfect for indoor use or in social situations. Navigate small and narrow spaces or corridors easily with its small frame. Although compact, the looks are grand!

POWERED WHEELCHAIR WITH ENDLESS POSSIBILITIES

With more than 20 colour options, a large number of options, and adjustments.... Forest 3 can be personalized to any driver and their preferences. Although every Forest is different, the quality and durability remain constant. Every user story is unique and we gladly adapt our chair to suit yours!



Forest 3

ADULT POWERED
WHEELCHAIR

C8
yellow

C53
orange

C85
green



We offer more than 32
different colours.
For more information,
please check our website.



SCAN the QR code and
enjoy the Forest 3 movie.

Compact, durable, and high performing. Forest 3 has all the best features for users looking for a compact powered wheelchair, suitable for in- and outdoor use. Combining strength, comfort, and adaptability, Forest 3 has a well-thought-out design that combines all of our expertise of more than 55 years.

The front-wheel axle is an undeniable asset to a powered wheelchair. It guarantees constant contact of the 4 wheels


with the ground whilst remaining stable and comfortable. The V-stab technology (optional) allows the driver to maintain a straight course when faced with rougher terrain or obstacles without constantly having to use compensating steering techniques.

For daily use, we all can agree that comfort is an important factor. With its preshaped backrest and individual adaptations, Forest 3 provides a personalized solution.

Features



*Electric-powered inclination
of the backrest and seat*

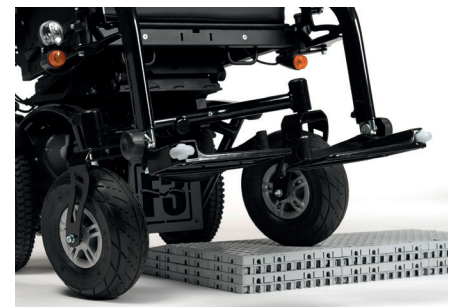

*Rotating axle to allow
constant ground contact
on all the wheels*



*Adjustable and foldable
armrests with memory
position*



OPTIONS Forest 3



B66MAL Guttered rotating armrest



Armrest preshaped PU



Armrest artificial leather upholstery



Armrest soft upholstery



B06 Removable and swing away footrests



BZ8 ART Articulating comfort legrests



BZ8E ART Powered articulating comfort legrests



L34 Pre-shaped backrest (soft)



L35 Pre-shaped seat (soft)



L16 + L35 Removable seat cushion with fasteners + board



TUB + L34 Tubular pre-shaped backrest with adjustable aluminium clamps



TUB + L44 Tubular backrest with adjustable tension




SE53 Swing away controller



S53A Swing away controller with lock positions



S53Z Swing away down controller for easy table access

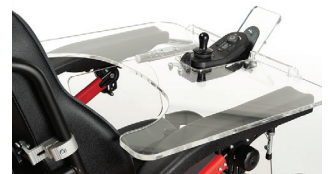

SE02 Controller integrated in the table



LINX LCD Electronic display



SE04 Chin controller



SE08 Controller at the back



SE09 Dual control accompanying person

Variations

Forest 3+

Compact, solid and efficient. The **FOREST 3+** is the ideal companion for users who are looking for a **compact wheelchair** for indoors and at the same time a high-performance model for outdoors. Strength, comfort and adaptability are the main features of the **FOREST 3+**.



Products designed to suit you

Our basic principle is: the wheelchair must be adapted to you, not the other way around. That is why our wheelchairs are evolutive and an answer to your specific needs.

All types of wheelchairs can be adapted, have numerous options and can be fully personalized so that there is a wheelchair for every person. This is the result of our continuous research and experience since 1957.

Technical information

■ Total width	595 / 610 / 630 / 650 mm
■ Total length	1185 / 1205 / 1224 mm
■ Total height	1036 mm
■ Seat height	430 - 480 mm
■ Seat width	400 / 420 / 450 / 480 / 500 mm
■ Seat depth	420 - 520 mm

Forest 3 SU

The Forest 3 SU is a power wheelchair with a stand-up possibility. Just like its brother, the Forest 3 but with additional functions. Complete with a safety belt at chest level and a padded knee guard for complete comfort and safety. The control unit allows the user to shift from the standing position to the sitting position to the 180° reclined position.



Quality guaranteed

Regarding quality, the **Forest 3** complies with the highest norms, laid down by certified test institutes in most European countries.

■ Batteries	2 x 12V / 80Ah - 220W or 350W
■ Speed*	6 km/h or 10 km/h
■ Driving range*	42 or 50 km
■ Weight	139 kg
■ Maximum user weight	130 kg
■ Additional information	 

* The driving range is calculated according to European standards. * Speed and driving range may deviate (ca. 15%). Permitted dimensional and weight tolerances: 15 mm, 1.5 kg and 1.5°. Subject to technical alterations. Images for illustration purposes only. Certain features may be optional.

VERMEIREN UK LTD

Suite 1&2 Marshall Business Centre - Faraday Road
HR4 9NS Hereford - United Kingdom
Company registration number 16075613

www.vermeiren.co.uk
info@vermeiren.co.uk



Contact your local dealer for more information.

VERMEIREN
we care for you