

# Club UL

## Ultra-low bed

**Ultra-low beds are a great solution for people with dementia, disorientation or in need of bed in closer proximity to the floor.** Club UL is a high-low bed that can be lowered to the floor, preventing falls and injuries. It helps you get in and out of bed easily, reduce anxiety and stress, and improve your independence and quality of life.



**VERMEIREN**  
we care for you

# Club Ultra Low (UL)

**A series of beds combining over 65+ years of expertise in home care solutions.** The Club UL from **VERMEIREN** is designed based on input from users, rental companies and health care professionals. A simple and safe four-section bed to provide the best comfort and care at home.

Adjust the Club UL bed electrically with their three Linak™ standard motors. Adjust the height according to the needs of either user or caregiver for care or transfers.

The Anti-Trendelenburg position is standard included to go the extra mile for the user and bring additional health benefits when needed.



The controller is easy to use for both carer and user but can be locked if needed for additional safety.



All 4 wheels can be locked for safety and security and allow for easy maneuvering.

## Patientlifts

**VERMEIREN** offers patientlifts that are compatible with our ultra-low beds. A patientlift is a device that helps to transfer patients from one place to another, such as from a bed to a wheelchair or from a wheelchair to a toilet. Club UL can be lowered to the floor level or raised to a comfortable height for the patient and the caregiver.



To ensure the best solution for you and your comfort, the Club UL bed has been designed with many configuration options. From the slatted base to the finishing covers.

**DISCOVER MORE Club UL features below.**



## Slatted base



Metal slats




Wooden slats

## Optional rastomat



Without the rastomat

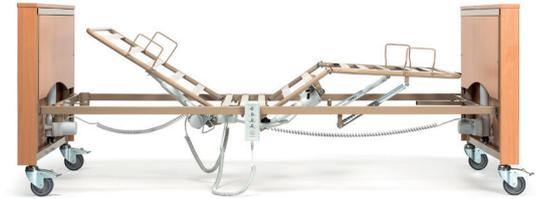

With the rastomat

## Side rails



Bed side rails in the upper position providing comfort and security



Adjust the side rails to the lower position for easy transfers or nursing

## Benefits of the ultra-low bed

It can prevent or reduce the impact of falls for patients who are at risk of falling out of their beds. For example, those with dementia, Parkinson's, epilepsy, or other conditions that cause disorientation, involuntary exits of the bed and general anxiety

- It can improve the comfort and dignity of the patients by allowing them to get in and out of the bed more easily and independently.
- It can enhance the safety and efficiency of the caregivers by reducing the physical strain and workload of lifting and transferring patients.




## Options

We've selected the best **Club UL** options for you. Discover all options here.



A monkey pole adds additional support to the independency of patients when lifting or repositioning themselves.

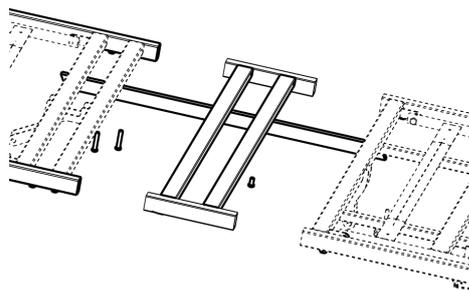

Club beds are created to suit everyone's needs! Easily add an additional length of 20 cm to the Club bed with our bed extension.



Serum holder



Transport kit

Configurations may differ depending on your region. Please consult our order forms for the most up-to-date information regarding all possibilities.

## Products designed to suit you

**VERMEIREN** always strives to come up with the best solution available to encompass your personal needs. As technology advances, so do we and as such we can offer you the most up to date products. This is the result of our continuous research and experience since 1957.

## Quality guaranteed

Regarding quality, **Club UL** complies with the highest norms, laid down by certified test institutes in most European countries.

## Technical data Club UL

■ Total width	1037 mm	■ Maximum angle for upper legs	20,7°
■ Total length	2175 mm	■ Maximum angle for lower legs	-16,5°
■ Total height	725 - 1135 mm	■ Max. backrest adjustment	69°
■ Height of mattress support (without mattress)	215 - 625 mm	■ Weight	73,5 kg
■ Mattress spec.	RG 35, thickness 120 mm, 2000x900 mm	■ Maximum user weight	165 kg
■ Anti-Trendelenburg	yes	■ Additional information	CE

Measurements in mm and kg. Deviations in sizes and weight (15 mm and 1,5 kg) are possible. Subject to technical alterations. Images for illustration purposes only. Certain features may be optional.

## VERMEIREN UK LTD

Suite 1&2 Marshall Business Centre - Faraday Road  
HR4 9NS Hereford - United Kingdom  
Company registration number 16075613

[www.vermeiren.co.uk](http://www.vermeiren.co.uk)  
[info@vermeiren.co.uk](mailto:info@vermeiren.co.uk)



For more information, contact your local dealer

**VERMEIREN**  
we care for you